Table Of Contents

Precautions	2
Contraindications	5
Warnings	5
DRF™ (Dynamic radio frequency) Technology – How does it work?	7
Intended Use	8
Possible Adverse Effects	9
Device Description	10
User interface description	11
User Instructions	12
Preliminary Test	12
Getting Started	13
Step 1: Prepare your skin	13
Step 2: Device Operation	13
Recommended Treatment	15
Post facial treatment care	17
Cleaning & Storage	18
Troubleshooting	18
Disposal	19
Customer Service	19
Device specification parameters	20
Labeling	21

Carefully read the information below before you begin your first treatment. If you still have any doubts or are unsure that the device is safe for your use, please consult with your physician.

Precautions

- Never use flammable liquids, such as alcohol (including perfumes, sanitizers, or other applications containing alcohol) or acetone to clean the skin in the hour before using the device, as it may cause burning of the skin.
- Do not open or repair the device, as it may expose you to dangerous electrical components.
- Do not place, store or use the device near bathtubs, showers, basins or

other vessels containing water, as this may cause severe electrocution.

- Do not leave the device unattended while plugged into the electricity supply, as it may lead to electrocution.
- Keep out of reach of children and pets.
- Not to be used by persons under 18 years of age.
- Children must not play with the appliance. Cleaning and maintenance by the user must not be carried out by children without supervision.
- Use the device only in accordance with these Instructions for Use and for its intended use.

- Do not use the device if any part of the device is missing, damaged or lost, if the device is not working properly, if you see or smell smoke, or if the device becomes damp or wet.
- Make sure the device is clean, dry and intact before each use. Cleaning instructions are detailed in the **Cleaning and Storage** section below.
- Do not use the device without first applying the Base Gel.
- Do not use the Base Gel after its expiration date.

Note: The Base Gel is for external use only – avoid direct contact with eyes. In the event that such contact does occur, wash the eyes immediately with water.

Contraindications

Do NOT use the device if you:

- Have an active implant, such as a pacemaker or insulin pumps, or if you have any known severe heart disorder, such as arrhythmias.
- Are pregnant or breastfeeding.
- Suffer from severe or uncontrolled endocrine or metabolic conditions, such as diabetes or uncontrolled thyroid disease.
- Have a damaged immune system or autoimmune disorder, such as AIDS or Lupus, or if you use drugs that influence the immune system.
- Suffer from blood clotting disorders, excessive bleeding or bruising, damaged wound healing or if you routinely take blood-thinning medication.
- Suffer from a history of deep vein thrombosis.
- Have a history of collagen disorders or suffer from keloid scarring (growth of extra scar tissue where the skin has healed after an injury).
- Have a general extreme illness, such as hepatitis, liver disease or epilepsy.
- Suffer from sensitive or damaged skin, such as in cases of frequent redness or rash, inflammations, breakouts (such as psoriasis), blisters, sores, infections, acne or extra dry skin.
- Have diminished or exaggerated perception of temperature change.
- Have taken Accutane in the past 6 months.
- Have a heat-induced disease or an active skin disease, such as herpes, eczema, rashes etc., in the treatment area.
- Have skin cancer, areas at risk of malignancy, or if you have received radiation therapy or chemotherapy within the past 3 months.

Warnings

Do NOT use the device:

- If you have any permanent implants in the designated treatment area, including but not limited to: silicon or injected chemical materials, metal plates or other metal materials.
- If you had plastic surgery, such as face lifting or eyelid surgery, or if you had aesthetic injections of fillers such as: Hyaluronic Acid, Collagen, fat or Botulinum Toxin, in the treatment area within the past 3 months.
- If you have had a professional chemical peel within the past 3 months with a Hydroxyl Acids (AHA), such as glycolic acid, or β Hydroxyl Acid (BHA), or if you had a professional facial abrasion such as: dermabrasion, facial resurfacing, deep chemical peeling (phenol, Trichloroacetic acid-TCA), and/or ablative laser treatments within the past 3 months.
- Over nostrils, lips of the mouth, ears, breast or genital area, or areas not recommended in this user manual.
- Over the thyroid since it may cause adverse events.
- Over the eyelids, eyes or sides of the eyes (next to your nose), since it may cause eye damage.
- If you have any metal items, such as jewelry or piercings, or tattoos or permanent make up on the designated treatment areas.
- Over sunburns, dark brown or black spots such as: freckles, birth marks, moles or blisters in the treatment area.
- If you are allergic or sensitive to chrome, nickel or any component of the Base Gel.
- The appliance may be used by persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge provided that they are under supervision or have been instructed in the safe use of the appliance and that they understand the hazards involved.

- Do not use the device over hairy areas of the body, as this may result in the hairs being burnt.

Note: If you have undergone any other aesthetic procedures in the designated treatment area, or if you have treated the skin with a prescribed drug within the past 3 months, or if you have any concerns regarding potential adverse effects, please consult your physician before using the device.

DRF™ (Dynamic radio frequency) Technology – How does it work?

This device brings a new approach to non-invasive skin rejuvenation featuring Sensica proprietary Dynamic RF technology, combining clinically proven RF (Radio Frequency), commonly used by doctors, with advanced facial massage.

Radio Frequency (RF) is a form of electromagnetic energy, intended to reduce wrinkles, tighten the skin, and renew skin elasticity. It is suitable for all skin tones, easy to use, and provides noticeable, long-lasting results.

RF technology has been used in doctors' clinics for several years and was found to be highly effective. The device has two electrodes that deliver the RF energy into the skin to contract collagen fibers and stimulate the natural production of collagen, which decreases as we age. If used correctly and according to the **Recommended Treatment**, the device will reduce the

appearance of wrinkles as well as tighten the skin and renew elasticity. Results may be visible after just one treatment.

Please make sure to carefully read this user guide before you start your treatment.

The user guide can also be found at www.sensica.com

Intended Use

The device is intended:

- For women and men who want to: reduce facial wrinkles, tighten their facial skin, and renew skin elasticity.
- For use on: face, cheeks, "smile" folds, "crow's feet", forehead, chin, jawlines, allowed areas of the neck, décolletage, and hands.
- For use on fine to moderate depth wrinkles.
- To be used with the Base Gel. This gel is especially designed to work with Sensica RF devices and ensures proper delivery of RF energy. This will provide the best possible results by improving contact between the device and the skin to improve skin conductivity.

Note: Sensilift is especially designed for your own personal treatment and is intended for a single user and home use only.

Possible Adverse Effects


When using the device according to the Instructions for Use, adverse effects are rare. However, every aesthetic procedure, including the usage of home-use devices, may cause some degree of adverse effects (as listed in the table below).

In case adverse effects appear, stop using the device and consult with your physician if necessary.

Mild adverse effects should be spontaneously relieved within a few hours.

Severity of Adverse Effect	Frequency of Adverse Effect	Adverse Effects
Mild	Common	Local pain Skin irritation Skin redness Skin sensitivity Allergies or sensitivity to the Base Gel
Moderate	Rare	Mild skin burns Pigment changes Excessive redness and swelling Urticaria, purpura or ecchymosis Hematoma
Severe	Extremely Rare	Severe burns Infection Bruising Blistering

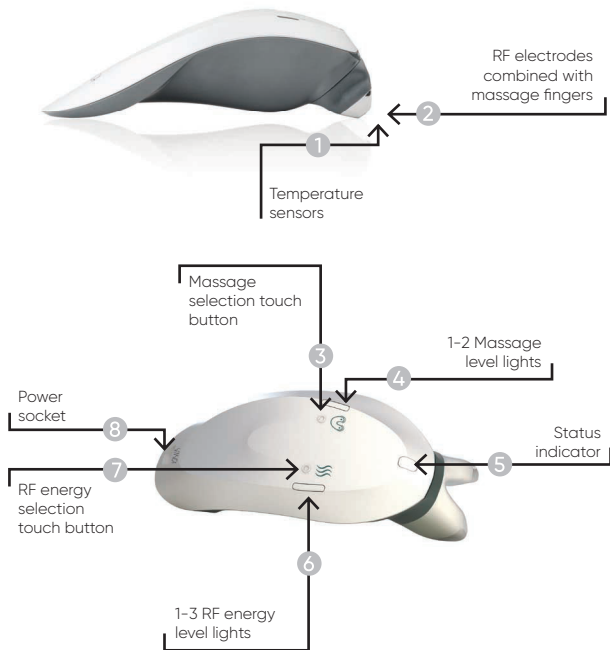Device Description

User interface description

1. **ON/OFF and energy selection button** – A gentle touch on the button (after the device is connected to the power supply) will turn the device ON at energy level 1, shown by 1 illuminated light in the energy level display. Once it reaches the highest level, another gentle touch will turn the device OFF. **The device will switch to standby mode when left unattended for 45 seconds.**
2. **Status indicators**

	Blinking yellow light	The device is in standby mode . In case of no light, make sure the power cord is plugged-in properly.
Status indicator	Steady green light	The device is in treatment mode . RF is activated, skin contact is detected, and the treatment is undergoing.
	Blinking alternating green and yellow lights	The device is in error mode . The device detected a malfunction and prevents the activation of the RF and massage functions.
RF energy level selection touch button 	1-3 energy levels lights	When the green light is on, it indicates that the RF energy is active. Touch the RF energy level selection button repeatedly until you reach the desired level: 1, 2, 3 or off.
Massage level selection touch button 	1-2 massage level lights	When the green light is on, it indicates that the massage function is active. Touch the massage level selection button repeatedly until you reach the desired level: 1, 2, or off.

User Instructions

The device is intended for use on: cheeks, nasolabial “smile” folds, “crow’s feet” (next to your eyes), forehead, chin, allowed areas of the neck, décolletage and hands.

What to Expect

During your treatment, it is normal to experience:

- Sensations of warmth created by the RF energy.
- Some mild red or pink color on the skin – during and just after each treatment.

Preliminary Test

The preliminary test is intended to determine your skin’s reaction to the treatment and help you find the optimal settings for each treatment area.

Perform a preliminary test on each area you intend to treat 48 hours prior to the first full treatment.

- Prepare your skin as described below in **Step 1: Prepare your skin.** Instructions for operating the device can be found below in **Step 2: Device operation.**

- Before using the device, apply an even layer of the Base Gel on the treatment area.
- Start your test with energy level 1 and massage level 1 and move the device in small circular motions.
- Increase the energy and massage levels according to your comfort. Stop the test if you do not feel comfortable or if you see adverse effects.
- Perform the preliminary test for a total time of about 3 minutes.
- Wait 48 hours and examine the tested area.
- If the skin in the tested area appears normal and no adverse effects are present, you can proceed with a full treatment at the highest comfortable energy level that you tested.

Note: You may select any RF energy level that you find comfortable. The RF energy level indicates the rate at which the temperature rises when in heating mode. A high RF energy level will enable

faster and more even heating. You can achieve maximum results with all RF energy levels, but, for lower levels, it will take more time to reach the effective temperature.

Getting Started

Step 1: Prepare your skin

1. Prior to use, make sure your skin is dry, clean and free of topical products such as: moisturizers, liquids or lint.
2. Remove jewelry from the treatment area as well as from your hands, prior to use.
3. Never use flammable liquids such as alcohol (including perfumes, sanitizers, or other applications containing alcohol) or acetone on the skin one hour before using the device as it may cause burning of the skin.
4. The device is suitable for men. However, male users should make sure their skin is clean-shaven before treatment.

Step 2: Device Operation

Note: A preliminary test should be conducted separately for each area you intend to treat before proceeding with a full treatment in that area. See **Preliminary Test** section above.

1. Make sure the device is clean, dry and intact prior to use.
2. Apply an even layer of Base Gel, included in your kit, on the area of your skin intended for treatment.
3. Set the Radio Frequency (RF) energy level and the massage speed to the level you selected in your preliminary test treatment. Gently touch the RF energy level selection button  repeatedly until you reach the desired level: 1, 2, or 3. Gently touch the massage level  button repeatedly until you reach the desired level: 1 or 2.
4. The activation of either RF or massage functions can be done independently or together. Both functions will be activated at level 1.

Select different levels based on the preliminary test.

5. Sensilift is equipped with 2 safety sensors which will keep the temperature in the effective and safe range during the entire treatment.

Note: The device stops the heating when desired temperature is achieved; heating is then re-activated when the temperature begins to drop.

6. Place the treatment fingers on the intended area, ensuring full contact with your skin during the entire treatment. Sensilift will indicate the proper skin contact by the green status light. If Sensilift is not properly placed, the RF energy will not be delivered.

Note: The device will automatically turn-off after 45 seconds of no skin contact.

7. If you do not feel sensations of warmth after approximately 10 seconds, try moving Sensilift on a smaller treatment area, and/or increasing the RF energy.
8. Move the device on the treated area in slow circular or back and forth motion while constantly and firmly contacting the skin.
9. See treatment areas, treatment time and treatment schedule in the **Recommended Treatment** section.
10. After completing the treatment, wipe off the remaining Base Gel from the treated area.
11. Cleaning instructions for before and after treatment are depicted in the **Cleaning & Storage** section.

Recommended Treatment

- Place the treatment fingers as indicated in the picture of the chosen treatment area.
- Continuously move the device as described for the chosen treatment area.

Treatment Area	Recommended Treatment Time	Recommended Device Movement
Forehead 	10 minutes	Slow up and down motion across the forehead
"Eleven" 	5 minutes	Small circular movements or up and down
"Crow's feet" 	5 minutes for each side	Continuous, slow 'side to side' circular motions
Cheeks 	5 minutes for each side	Slow circular motions across the cheek

"Smile folds"		5 minutes for each side	Slow up and down motions over the "smile folds"
Chin		5 minutes	Side to side
Jawline		5 minutes for each side	Slow up and down motions along the jaw
Neck		5 minutes for each side	Slow up and down motions across the neck
Decolletage		5 minutes for each 4cm ² sub-area	Slow circular motions across decolletage
Hands		5 minutes for each 4cm ² sub-area	Slow circular motions across hands

Warning: avoid contact with the eyes, eye cavities, nostrils and oral cavity, as treatment in these areas may cause adverse effects.

Warning: avoid treatment in the lower third of the neck since this area is more sensitive and avoid treating the thyroid area (marked Δ).

- For long-lasting results, use the device once a week on each desired area for a period of 8 weeks.
 - Response to the treatment differs between individuals and depends on environmental and biological factors, heredity and aging.
 - To sustain optimal long-term results, maintenance treatments may be performed as needed every 4–8 weeks.
- We recommend you also use Sensica's anti-aging cream to increase treatment efficacy and ensure long-lasting results (available for purchase at www.sensica.com).
 - Apply a high-factor sunscreen (at least 30 SPF) and protect the treated area from sunlight for two weeks following treatment.
 - Following the treatment, avoid any usage of hot water on your face and discontinue usage of abrasive or harsh products for two days.
 - Makeup may be applied immediately after treatment.

Post facial treatment care

Cleaning & Storage

- Make sure to turn off and unplug the device before cleaning.
- It is recommended to clean the device, especially the treatment fingers, temperature sensors, and silicone cover, of residual Base Gel with a soft damp cloth or wipes.

Do not rinse the device under running water, as it may cause device malfunction.

- Keep the device dry. Never place, store or use the device near: bathtubs, showers, basins or other vessels containing water as this may cause severe electrocution.

Troubleshooting

If you experience any problems with your device, please refer to the information below to attempt to resolve it:

- If the device does not turn on:
 1. Make sure that the power cord is plugged properly into the electrical outlet and into the device power socket.
 2. Make sure the power outlet has power.
- If you cannot feel heat:
 1. Verify full contact with the skin (indicated by a green light).
 2. Verify that you applied enough Base Gel.
 3. Remember that the device reacts to your skin and may temporarily reduce heat for safety reasons.

See status indicator table in the **User Interface Description** section.

- In case of dirt or debris on the treatment fingers, you should clean them in accordance with the **Cleaning & Storage** section. If you can't get the area clean, contact our customer service.
- In case of a general error, reset the device by unplugging the power cord from the electrical outlet on the wall, wait 10 seconds and then reconnect it.

If you still have problems with your device, please contact our customer service.

Disposal

For environmental reasons, do not dispose of the device in the household waste. Dispose it at a suitable local collection or recycling point in accordance with EC Directive–WEEE (Waste Electrical and Electronic Equipment). If you have any questions, please contact the local authorities responsible for waste disposal.

Customer Service

If your device has broken/ stopped working or if you need any further assistance/ information, please contact us. You may also contact us to order additional Base Gel or you can purchase it online at [sensica.com](https://www.sensica.com).

For customer service, please see contact details on inner back cover.

Device specification parameters

The RF handheld device is powered by bipolar RF at 10 ± 0.5 Watt

Frequency

1 ± 0.05 MHz

Selectable output power level

3 levels

Mechanical massage

Implemented by moving RF electrodes.

Speed

50–100 RPM

Temperature sensors

RF heating is automatically deactivated when skin temperature reaches $40^{\circ}\text{C}/104^{\circ}\text{F}$. The heating is reactivated when temperature decreases.

Treatment area

$\sim 2.2\text{cm}^2 / 0.34\text{in}^2$

External power supply

Input 100–240VAC, 50–60Hz,
Output 9VDC, 1.5A

Device size and weight

140g, 173mm x 56.8mm x 53.1mm
4.93oz, 6.8in x 2.23in x 2.09in

Storage Instructions

Storage and transport conditions

Temperature	-25°C to 70°C (-13°F to 158°F)
-------------	---

Humidity	10% to 90% RH at $55^{\circ}\text{C}/131^{\circ}\text{F}$
----------	--

Pressure	Non-condensing
----------	----------------

Operational conditions

Temperature	5°C to 40°C (41°F to 104°F)
-------------	--

Humidity	20% to 80% RH
----------	---------------

Max operating altitude	3,000 m.
------------------------	----------

Labeling

Model: STG100

Power: 9V, 1.5A, 100-240VAC,
50-60Hz

Double insulation

Consult Instructions for use

Waste Electronic and
Electrical Equipment

Please read all warnings
and safety precautions
prior to use

Protect from moisture.

This device must not be
used near to water or in
water itself (e.g., in a wash
basin, in the shower or
bathtub) – risk of electric
shock!

Our products and
packaging are recyclable,
do not throw them away!

Use the device only
with the original
power supply. For
your safety, never use
DC adapters which
are not approved by
the manufacturer.

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Warranty Card

Serial Number

Congratulations on your Sensica purchase!

We have every faith in the quality of each Sensica device - we know they meet the highest standards, because we subject them to rigorous quality and performance testing. In the unlikely event that you do find a defect in a Sensica device, it may be covered by this warranty.

Two Year Limited Warranty

This warranty covers only Sensica devices. Other Sensica products, such as pads or creams, are defined as consumables and therefore excluded from warranty coverage.

This limited warranty covers manufacturer defect related problems in your Sensica device, if you purchased it directly from www.sensica.com or from an authorized partner. Warranty approvals will only be granted

where an original receipt is provided. This warranty is limited to devices proven to be defective.

In order to obtain warranty service on your Sensica device, please go to our website or call our customer service center for assistance (you will find the number at www.sensica.com).

Sensica requires proof of purchase. Please provide an invoice or purchase receipt showing the date of purchase of the device

Enter your serial number on our website to activate your warranty.

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