## Table Of Contents

Precautions	2
Contraindications	5
Warnings	6
Introducing LIPOTHERM RF™ technology - How it works	7
Intended use	7
Possible Adverse Effects	8
Device description	9
User interface description	10
User instructions	11
Preparing the treatment area	11
Preliminary Test	12
Device Operation	13
Recommended treatment	14
Treatment areas	15
Post-treatment care	17
Troubleshooting	18
Device specification parameters	19
Labeling	20

## Precautions

- Never use flammable liquids, such as alcohol (including perfumes, sanitizers, or other applications containing alcohol) or acetone to clean the skin in the hour before using Sensifirm as it may cause burning of the skin.
- Do not open or repair the device, as it may expose you to dangerous electrical components.
- Do not place, store or use Sensifirm near bathtubs, showers, basins or other vessels containing water, as this may cause severe electrocution.
- Do not leave Sensifirm unattended while plugged into the electricity supply, as it may lead to electrocution.

- Keep out of reach of children and pets.
- Not to be used by persons under 18 years of age.
- Children must not play with the appliance. Cleaning and maintenance by the user must not be carried out by children without supervision.
- Use the Sensifirm only in accordance with these Instructions for Use and for the intended use.
- Do not use the device if any part of the device is missing, damaged or loose, if the device is not working properly, if you see or smell smoke, or if the device becomes damp or wet.
- Make sure the device is clean, dry and intact before each use. Cleaning

instructions are detailed in the Cleaning section below.

- Do not use the Sensifirm device without first applying the Base Gel.
- Do not use the Base Gel after its expiration date.

Note: The Base Gel is for external use only - avoid direct contact with eyes. In the event that such contact does occur, wash the eyes immediately with water.

## Who should NOT use the device?

Carefully read the information below before you begin a treatment with the device. If you still have any doubts or are unsure that the device is safe for your use, please consult with your physician.

## Contraindications

Do NOT use the device if you:

- have an active implant, such as a pacemaker or insulin pump, or if you have any known severe heart disorder, such as arrhythmias.
- are pregnant or breastfeeding.
- suffer from severe or uncontrolled endocrine or metabolic conditions, such as diabetes or uncontrolled thyroid disease.
- have a damaged immune system or autoimmune disorder, such as AIDS or Lupus, or if you use drugs that influence the immune system.
- suffer from blood clotting disorders, excessive bleeding or bruising, damaged wound healing or if you routinely take blood-thinning medication.
- suffer from a history of deep vein thrombosis.
- have a history of collagen disorders or suffer from keloid scarring (growth of extra scar tissue where the skin has healed after an injury).
- have a general extreme illness, such as hepatitis, liver disease or epilepsy.
- suffer from sensitive or damaged skin, such as in cases of frequent redness or rash, inflammations, breakouts (such as psoriasis), blisters, sores, infections, acne or extra dry skin.
- have diminished or exaggerated perception of temperature change.
- have taken Accutane in the past 6 months.
- have a heat-induced disease or an active skin disease, such as herpes, eczema, rashes etc. in the treatment area.

- have skin cancer, areas at risk of malignancy, or if you have received radiation therapy or chemotherapy within the past 3 months.

## Warnings

Do NOT use the device:

- If you have any permanent implants in the designated treatment area, including but not limited to: silicon or injected chemical materials, metal plates or other metal materials.
- If in the past 3 months you have had surgery, including plastic surgery, or you have undergone any aesthetic procedure, such as skin resurfacing, liposuction, deep peeling, cryolipolysis (fat freezing), or have had aesthetic injections of fillers, such as hyaluronic acid, collagen, fat or Botox, in the treatment area.
- If you have any metal items, such as jewelry or piercings, or tattoos or permanent make up on the designated treatment areas.
- If you have sunburn, or dark brown or black spots such

as freckles, birth marks, moles or blisters in the treatment area.

- If you are allergic or sensitive to chrome or nickel or any component of the Base Gel.
- The appliance may be used by persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge provided that they are under supervision or have been instructed in the safe use of the appliance and that they understand the hazards involved.
- Do not use Sensifirm over breasts, genitals or areas not recommended in this Instructions for Use.
- Do not use Sensifirm over hairy areas of the body, as this may result in the hairs being burnt.

---

**Note:** If you have undergone any other aesthetic procedure in the designated treatment area, or if you have treated the skin with a prescribed drug within the past 3 months, or if you have any concerns regarding potential side effects, please consult your physician before using Sensifirm.

---

## Introducing LIPOTHERM RF™ technology - How it works

Radio Frequency (RF) is a form of electromagnetic energy that has been used for a long time in aesthetic medicine. The spacing between Sensifirm's four RF electrodes enables great treatment volume and allows multiple RF waves to reach deeply into the skin layers, heating the dermis, hypodermis (inner skin layers) and subcutaneous fat layers, without damaging the epidermis (outer skin layer). This deep thermal effect in turn activates several biological mechanisms, among them remodelling, and later producing, collagen and elastin; decreasing the volume of adipocytes' (fat cells); and improving microcirculation and lymphatic drainage. Together, these effects lead to firming and smoothing of the skin, improvement in cellulite appearance and circumferential reduction, for a tighter body with the best possible visible, lasting results.

Please make sure to carefully read these Instructions for Use before you start your Sensifirm treatment. This guide, as well as an instructional video, can also be found at [www.sensica.com](http://www.sensica.com)

---

## Intended use

- By women who want to improve cellulite appearance.
- By men and women who want to tighten and reduce the circumference of their body.
- Sensifirm can be used on buttocks, thighs, abdomen, waist and arms.

---

**Note:** This device is especially designed for your own personal treatment and is intended for a single user and home use only.


---

## Possible Adverse Effects

When using this device according to the Instructions for Use, adverse effects are rare. However, every aesthetic procedure may cause some degree of adverse effects (as listed in the table below). In case adverse effects appear, stop using the device, and consult with your physician if necessary. Mild adverse effects should be spontaneously relieved within a few hours.

<b>Severity of Adverse Effect</b>	<b>Frequency of Adverse effect</b>	<b>Adverse effects</b>
<b>Mild</b>	<b>Common</b>	Local pain Skin irritation Skin redness Skin sensitivity Allergies or sensitivity to the Base Gel
<b>Moderate</b>	<b>Rare</b>	Mild skin burns Pigment changes Excessive redness, rash, swelling Subcutaneous bleeding
<b>Severe</b>	<b>Extremely rare</b>	Severe burns Infection Bruising Blistering

## Device description
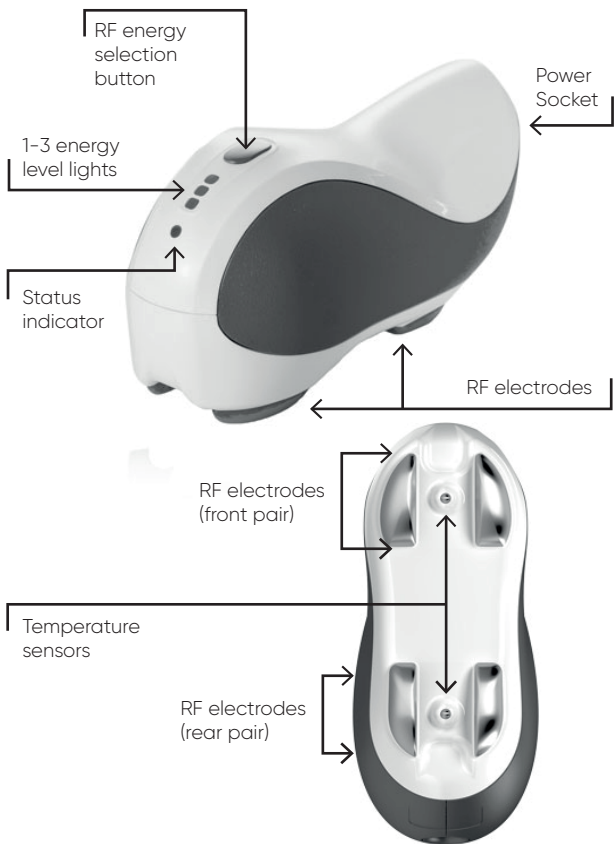
## User interface description

Status Indicator	Transient orange light (for ~5 sec) followed by transient green light (for ~1 sec)	Sensifirm is in <b>standby mode</b> and ready to work.
	Steady green light	Sensifirm is in <b>treatment mode</b> . RF is activated, skin contact is detected and the treatment is underway*.
	Alternating red and orange lights	Sensifirm is in <b>error mode</b> . The device has detected a malfunction and has blocked the activation of the RF function.
RF energy level selection button	1-3 energy levels lights	The green light indicates which RF energy level has been activated. Press the RF energy level selection button repeatedly until you reach the desired level: 1, 2, or 3.

\* **Note:** RF energy is delivered only when proper contact with your skin is detected. When treatment is underway, the status indicator will have a steady green light and you will feel vibration.

\* Upon contacting the treatment area, the massage function is automatically turned on.

## User instructions

Sensifirm is intended for use on buttocks, thighs, abdomen, waist and arms.

### What to expect when using Sensifirm

During your treatment, it is normal to experience:

- Sensations of warmth created by the RF energy.
- Some mild red or pink colour on the skin, during and just after the treatment session.

### Preparing the treatment area

1. area clean and free of topical products, such as moisturizers and lotions. Caution: Never use flammable liquids, such as alcohol (including perfumes, sanitizers, or other applications containing alcohol or acetone) to clean the skin in the hour before using Sensifirm, as it may cause burning of the skin.

2. Remove all jewelry from the treatment area and from your hands.
3. Remove long and/or thick hair from the treatment area (mostly applicable for male users).
4. Perform a preliminary test on a small patch of skin in the intended treatment area, 48 hours before the first use of Sensifirm. In the preliminary test, you can determine the energy level settings that best suit you (see Preliminary test section below).

---

#### Important:

- Check the expiration date of the Base Gel.
  - Before beginning the treatment, make sure that the electrodes are clean, dry and intact.
-

## Preliminary Test

The preliminary test is intended to determine your skin's reaction to treatment with the device and help you set the optimal setting for each treatment area.

**Perform a preliminary test on each area you intend to treat 48 hours before the first full treatment.**

1. Prepare the treatment area as detailed above.
2. Apply a thin, even layer of Base Gel, included in your kit, on the treatment area.
3. Plug in the Sensifirm device and select RF energy level 3.
4. Place the Sensifirm device on the desired treatment area and move it in a circular motion, for 3 minutes.
5. The treatment will create a warm sensation on your skin. However, if you feel your skin is heating up too fast, reduce the RF energy level, adjusting it to suit your own comfort.
6. Stop the test immediately if you do not feel comfortable, or if you see adverse effects such as severe redness or blistering.
7. Wait 48 hours and examine the tested area.
8. If the skin in the test area seems normal (no adverse effects) you can proceed with a full treatment at the highest comfortable energy level.

---

**Note:** We recommend performing all treatments at energy level 3 (highest). However, you can achieve maximum results at all RF energy levels: at lower levels, it will just take longer to reach the effective temperature. The RF energy level indicates the rate at which the temperature rises when in heating mode. A high RF energy level will enable you faster and more even heating.

---

## Device Operation

---

**Note:** A preliminary test should be conducted before proceeding to a full treatment. (See **Preliminary Test** section above).

---

1. Prepare the treatment area as detailed in “preparing the treatment area” section above.
2. Apply an even layer of the Base Gel included in your kit, on the area that is intended for treatment.
3. Set the RF energy level to the level you found most comfortable when you carried out the preliminary test. Press the RF energy level selection button repeatedly until you reach the desired level: 1, 2, or 3.
4. Place the device on the intended treatment area, ensuring full contact of the electrodes with your skin throughout the entire treatment. If the device is not properly placed, the RF energy and massage will not be delivered. When the device is working, you will feel and hear a gentle vibration, and the indicator will be lit up in green.
5. Move the device over the treatment area in slow circular motions, while maintaining constant, firm contact with the skin.
6. If you do not feel a sensation of warmth after approximately 30 seconds, try moving the Sensifirm over a smaller treatment area, and/or increasing the RF energy.
7. At the end of each treatment session, turn off the device by pressing on the energy selection button. An orange indicator light will light up while the device is shutting down. After the light has turned off, unplug the device from the electrical outlet.

8. After completing the treatment, wipe off any remaining Base Gel from the treated area.
9. Clean the device with a soft damp cloth or paper towel, according to the cleaning instructions in the **Cleaning** section below.

## Recommended treatment

For long-lasting results, use Sensifirm once a week on each treatment area, for a period of 8 weeks.

To sustain optimal results, you can repeat the Sensifirm treatment protocol as needed.

You can use the Sensifirm device on the buttocks, thighs, abdomen, waist or arms. Treatment session depends on area size.

1. 20 minutes for a body area the size of 15cmX30cm/ 6inX12in (this is about the size of your two palms).
2. 10 minutes for a body area the size of 8cmX15cm/ 3inX6in (this is about the size of one palm).


## Treatment areas

Do not treat full areas with energy levels you have not tested before (as detailed in the Preliminary test section above).

### Buttocks (area 1)

---

- Place the electrodes on the treatment area, as indicated in the diagram.
- Move the Sensifirm in continuous circular motions around the area, as indicated by the arrows in the diagram.
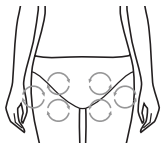
### Thighs (area 2)

---


- Place the electrodes on the treatment area, as indicated in the diagram.
- Move the Sensifirm in continuous circular motions around the area, as indicated by the arrows in the diagram.

### Abdomen (area 3)

---

- Place the electrodes on the treatment area, as indicated in the diagram.
- Move the Sensifirm in continuous circular motions around the area, as indicated by the arrows in the diagram.
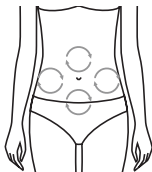
### Waist (area 4)

---


- Place the electrodes on the treatment area, as indicated in the diagram.
- Move the Sensifirm in continuous circular motions around the area, as indicated by the arrows in the diagram. If this is uncomfortable, switch to linear motions.

### Arms (area 5)

---

- Place the electrodes on the treatment area, as indicated in the diagram.
- In this area you can use just 2 electrodes.
- Move the Sensifirm in continuous circular motions around the area, as indicated by the arrows in the diagram. If this is uncomfortable, switch to linear motions.
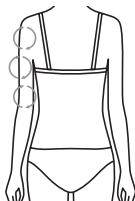
## Post-treatment care

- After completing the treatment, wipe off the remaining Base Gel from the treated area.
- Do not expose the treatment area to direct sun light without applying a high-factor sunscreen (at least 30 SPF). Uncontrolled sun exposure might cause your skin to be sensitive.
- For 2 days after treatment, avoid the use of extremely hot water and do not use abrasive or harsh products.
- Since skin cellulite formation is an ongoing process, you may repeat the Sensifirm treatment according to your needs.

## Cleaning

1. Make sure to turn off and unplug the device before cleaning.
2. It is recommended to use a soft damp cloth, paper towel or pre-moistened wipe to clean any residual Base Gel from the Sensifirm. Make sure the entire underside of the device, including the electrodes area, temperature sensors and plastic cover, is cleaned thoroughly.

Do not rinse the Sensifirm device under running water, as this may cause device malfunction.

---

## Troubleshooting

If you experience any problems with your Sensifirm, please refer to the information below:

<p><b>If the device does not turn on</b></p>	<ul style="list-style-type: none"> <li>• Make sure that the power cord is plugged in properly, both to the electrical outlet and into the device socket.</li> <li>• Make sure that the power outlet has power.</li> </ul>
<p><b>If Sensifirm does not deliver RF signal energy and massage (i.e. you cannot feel the warm sensation or feel the gentle vibration)</b></p>	<ul style="list-style-type: none"> <li>• Make sure that there is full contact of the electrodes with the skin.</li> <li>• Gently press the device onto your skin.</li> <li>• Make sure that there is enough Base Gel on the treatment area.</li> <li>• The RF signal will be activated only when at least 2 electrodes are in full contact with the skin. (See indicator details set out in the User Interface description section above.)</li> </ul>
<p><b>In case of general error</b></p>	<ul style="list-style-type: none"> <li>• Reset the device – unplug the power cord from the electrical outlet on the wall, wait 10 seconds and then reconnect it.</li> <li>• If dirt or debris get stuck on the electrodes, you should clean them in accordance with the <b>Cleaning</b> section above. If you cannot get the area clean, contact our customer service.</li> <li>• If you still have problems with your device, please contact our customer service.</li> </ul>

---

## Customer Service

For customer service, please see contact details on inner back cover.

---

## Device specification parameters

The RF handheld device is powered by multipolar RF at  $10\pm 0.5$  Watt

### Frequency

$1\pm 0.05$  MHz

### Selectable output power level

3 levels

### Temperature sensors

RF heating is automatically deactivated when skin temperature reaches  $42^{\circ}\text{C}$  ( $106.7^{\circ}\text{F}$ ). The heating is reactivated when temperature decreases.

### Treatment area

$\sim 10\text{--}20\text{cm}^2$  ( $\sim 3.25\text{--}7.5\text{in}^2$ )

### External power supply

Input 100–240VAC,  
50–60Hz, 0.5A  
Output 12VDC, 1.8A

### Device size (LxHxW)

146mm X 74mm X 63mm  
5.75in X 2.91in X 2.48in

### Weight

192g (6.77oz)

## Storage Instructions

---

### Storage and Transport Conditions

---

Temperature	$-10^{\circ}\text{C}$ to $70^{\circ}\text{C}$ ( $14^{\circ}\text{F}$ to $158^{\circ}\text{F}$ )
-------------	--

---

Humidity	10% to 90% RH at $55^{\circ}\text{C}$ / $131^{\circ}\text{F}$
----------	--

---

Pressure	Non-condensing
----------	----------------

---

### Operational Conditions

---

Temperature	$10^{\circ}\text{C}$ to $40^{\circ}\text{C}$ ( $50^{\circ}\text{F}$ to $104^{\circ}\text{F}$ )
-------------	---

---

Humidity	30% to 70% RH
----------	---------------

---

Max Operating Altitude	3,000 m.
------------------------	----------

---

Use Sensifirm only with the original power supply. For your safety, never use DC adapters which are not approved by the manufacturer.

---

## Labeling

Model: CEL100

Double insulation

Consult Instructions for use

Waste Electronic and Electrical  
Equipment

Please read all warnings and safety  
precautions prior to use

Do not rinse Sensifirm under running  
water, as it may cause device  
malfunction. Keep dry

Our products and packaging are  
recyclable, do not throw them away!

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









## Treatment log

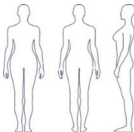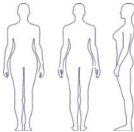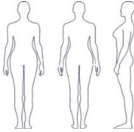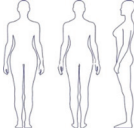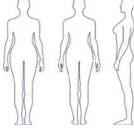Treatment date	Treated area	Comments
		
		
		
		
		

## Treatment log

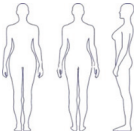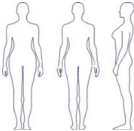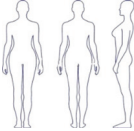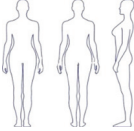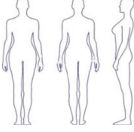Treatment date	Treated area	Comments
		
		
		
		
		

# Warranty Card

Serial Number

Congratulations on your Sensica purchase!

We have every faith in the quality of each Sensica device - we know they meet the highest standards, because we subject them to rigorous quality and performance testing. In the unlikely event that you do find a defect in a Sensica device, it may be covered by this warranty.

## Two Year Limited Warranty

This warranty covers only Sensica devices. Other Sensica products, such as pads or creams, are defined as consumables and therefore excluded from warranty coverage.

This limited warranty covers manufacturer defect related problems in your Sensica device, if you purchased it directly from [www.sensica.com](http://www.sensica.com) or from an authorized partner. Warranty approvals will only be granted

where an original receipt is provided. This warranty is limited to devices proven to be defective.

In order to obtain warranty service on your Sensica device, please go to our website or call our customer service center for assistance (you will find the number at [www.sensica.com](http://www.sensica.com)).

Sensica requires proof of purchase. Please provide an invoice or purchase receipt showing the date of purchase of the device

Enter your serial number on our website to activate your warranty.

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